

MK-3CN HEMATOLOGY CONTROL+CRP

2.5mL

**NORMAL
ASSAY SHEET**

LOT No.	MK3C22YN
EXP. DATE	2023-02-05
OPEN VIAL STABILITY	14DAYS

Manual Input

For manual input, enter the following on the new lot registration screen.

(This manual input is not required when entering the barcode.)

Lot	22Y
Level	Normal
Control	MK-3C

MEK-1303 Barcode Entry

Please read all three barcodes below.

If the barcode is difficult to scan, place this assay sheet on a flat surface, keep the barcode reader and barcode vertical, and scan within a distance of 2cm.


1



2



3

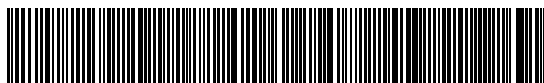

PARAMETER	UNIT	MEK-1303	RANGE
		ASSAY VALUES	
WBC	10 ³ /μL	7.65	± 0.80
RBC	10 ⁶ /μL	4.71	± 0.20
HGB	g/dL	13.82	± 0.50
HCT	%	45.6	± 3.0
MCV	fL	96.9	± 7.0
MCH	pg	29.3	± 2.4
MCHC	g/dL	30.3	± 3.0
PLT	10 ³ /μL	300.9	± 50.0
MPV	fL	9.1	± 2.0
RDW-CV	%	16.9	± 3.0
LY%	%	34.38	± 10.00
MO%	%	2.44	± 5.00
GR%	%	63.18	± 10.00
LY	10 ³ /μL	2.63	± 1.00
MO	10 ³ /μL	0.19	± 0.40
GR	10 ³ /μL	4.83	± 1.00
DILUENT		ISOTONAC·3 ISOTONAC·4	
HEMOLYSING REAGENT		HEMOLYNAC·310	
CRP ※1	Precision	10.0% or less	
	Reference range	17.70-32.90 mg/L	
	Measurement cell	Celltac chemi CRP 4N	

※1 Accuracy cannot be monitored by the CRP value of MK-3CN.
Use the CRP value only for precision (reproducibility) monitoring.

* The MEK-1303 can automatically input the assay values and ranges by reading all the barcodes 1, 2, and 3 on the left with the attached handy barcode reader. Barcodes 1, 2 and 3 can be read in any order.

* Please note that the assay ranges are reviewed lot by lot and the ranges are subject to change.

 NIHON KOHDEN CORPORATION 1-31-4 NISHIOCHIAI, SHINJUKU-KU TOKYO 161-8560 JAPAN
--



0694-901766